

ClariSure GK Polyethersulfone Membrane Cartridge Filter with Microglassfiber prefilter



mdi polyethersulfone (PES) membrane cartridge filters *ClariSure GK* are disposable filtration devices for clarification of heavily particle laden solutions or colloidal suspensions. There are especially suited for serum and plasma applications. These cartridges contain an upstream microglassfiber layer and a polyethersulfone membrane layer on the downstream.

The upstream microglassfiber layer retains large size particulate contamination and colloidal particles to protect premature clogging of downstream side PES membrane layer. Thus it gives high through-puts resulting in better economics.

mdi *ClariSure GK* filters are validated for key performance parameters such as chemical compatibility, extractables, heat stability and flow rates. These are available in a variety of sizes to suit specific microfiltration needs in critical and specialized process applications.

Special Features

- Large filtration area
- High throughputs
- Long service life
- Multiple autoclavable
- Heat sealed, no glues adhesives
- Low protein binding

Applications

- Prefiltration of parenterals
- Prefiltration of ophthalmic solutions
- Prefiltration of proteinaceous liquid where minimum protein loss is desired, such as sera, culture soups and recombinant proteins, antibodies etc.
- Prefiltration of media, buffers etc for particulate contamination reduction

Complies with USFDA 21 CFR 210.3(b)(6)


Meets and Exceeds USFDA 21 CFR 177.1520

Specifications

Construction				
Pore Size	0.1 µm, 0.2 µm and 0.5µm			
Filter Media	Hydrophilic PES			
Support Layers	Polyester			
Body and Core	Polypropylene			
Size				
Size	5"	10"	20"	30"
Effective Filtration Area (Nominal)	2500cm ²	5000cm ²	10000cm ²	15000cm ²
Operational				
Max. Operating Temperature	80 °C @ < 30 psi (2 Kg/cm ²)			
Max. Differential Pressure	50 psi (3.5 Kg/cm ²) @ 25 °C			
Sterilization	Autoclavable/Steam Sterilized at 121°C for 30minute, 25 Cycles.			
Assurance				
Toxicity	Passes Bioreactivity test, In Vivo, as per USP <88> for Class VI plastics			
Non Fiber Releasing	Passes test as per USP and comply with USFDA 21 CFR Part 210.3(b)(6) for fiber release			
Extractables with WFI	Passes test as per USP			
Oxidizable Substances	Within limits as specified in USP			
Particle Shedding	Passes USP test for particulates in injectables			
Indirect Food Additive	All Polypropylene components meet the FDA Indirect Food Additive requirements cited in 21 CFR 177.1520			
Good Manufacturing Practice	These products are manufactured in a facility which adheres to Good Manufacturing Practices.			
Quality Management System	ISO-9001:2008 Certified			
USFDA	DMF No. 015554			

Water Flow Rates

0.5 µm ClariSure GK Cartridge Filters

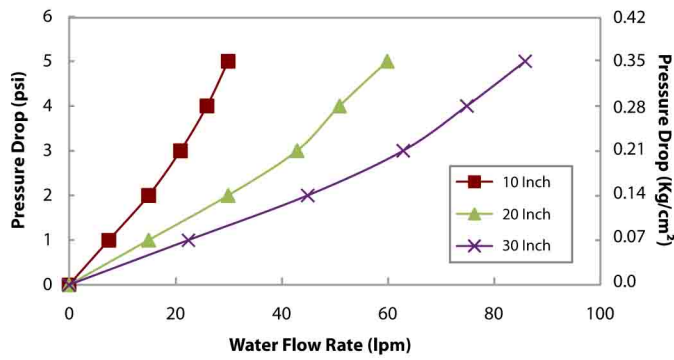

Ordering Information

Type		Size		Pore Size		Adaptor		Elastomer		Sterility		Pack Size	
	Code		Code		Code		Code		Code		Code		Code
Clarisure GK (Polyethersulfone Membrane Cartridge with Microglassfiber Upstream Layer)	CGKX	5"	53	0.1 µm	36	7P	A0	Silicone	SS	Non Sterile	1	1	01
		10"	54	0.2 µm	01	7P without Fin	A1	EPDM	SE				
		20"	55	0.5 µm	04	'0'	D0	Viton	SV				
		30"	56					FEP Encapsulated Viton	FV**				
Example													
CGKX		54		04		A0		SS		1		01	

* Size 5" are available in code A0 (7P) and A1 (7P without Fin) only

**FV is available in Adapter Code A0 (7P) and A1 (7P without fin) only



Gesfilter, S.L.

Distribuidor Oficial España

Avenida Cerdanyola, 92-94 Local B4

08173 - Sant Cugat del Vallès - Barcelona

Tel: 935921554 - Email: comercial@gesfilter.com

www.gesfilter.com

Advanced Microdevices Pvt. Ltd.

20-21, Industrial Area, Ambala Cantt-133 006, INDIA

Tel: +91-171-2699290, 2699471 Fax: +91-171-2699221

E-mail: info@mdimembrane.com Internet: www.mdimembrane.com